



AN UPTOWN COMBO AT A DOWNTOWN PRICE.

The Sonicator® Plus 930 Combo

Introducing a hot new combo that performs muscle stimulation *and* therapeutic ultrasound. The exciting new Sonicator Plus 930 features a 1- and 3-MHz, 5 cm² applicator for therapeutic ultrasound, plus the three most popular stim waveforms—Interferential, Premodulated and Medium Frequency (Russian).

In one compact package, the Sonicator Plus 930 gives you:


- Independent control of up to three treatment protocols
- Two intensity controls and two intensity displays, enabling the operator to monitor all aspects of a treatment simultaneously
- An easy-to-use membrane panel designed for quick set-up and simple operation
- A convenient recessed cradle that holds the ultrasound applicator between treatments

So contact Mettler Electronics today at 800-854-9305 to learn more about the value-priced Sonicator Plus 930. In one console, you get muscle stimulation *plus* therapeutic ultrasound. And all that jazz.

METTTLER ELECTRONICS®

©2001, Mettler Electronics Corp. All rights reserved.

THE SONICATOR® PLUS 930 COMBO

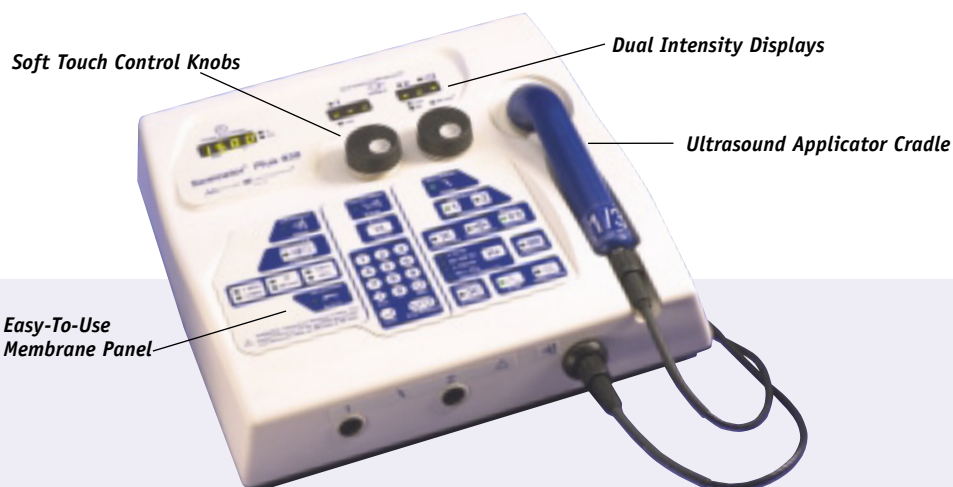

The model 93 or 94 cart offers a sturdy treatment platform for heavy clinical use.



The economical model 73 cart includes 3 shelves to accommodate Mettler units and accessories.

THE RIGHT COMBINATION OF PERFORMANCE AND VALUE.

The value-priced Sonicator Plus 930 combines a wealth of important features for the busy practice. Therapeutic ultrasound is provided by a 1- and 3-MHz, 5 cm² applicator with unique direct crystal-to-patient contact. And muscle stimulation is delivered in a choice of three waveforms—Interferential, Premodulated, and Medium Frequency (Russian stimulation).

The Sonicator Plus 930 is also remarkably easy to set up and use. To start a treatment, just select treatment parameters with as few as three keystrokes. Then adjust the intensity with a twist of the dial.

Impressive performance capabilities, outstanding simplicity of design—that's the Sonicator Plus 930. Dollar for dollar, no one gives you more than Mettler.

General Specifications:

Input:	90-240 VAC, 50-60 Hz
ETL and C-ETL Listed:	Model ME 930 (9801427)
Weight:	8.7 pounds
Dimensions:	6 in (H) x 12 in (W) x 12 in (D)
Warranty:	2 year limited

Ultrasound Specifications:

Effective Radiating Area:	5 cm ²
Frequency:	1.0 MHz 3.2 MHz
Modes:	Continuous Pulsed, 20%
Maximum Output Power:	11 W
Maximum Intensity:	2.2 W/cm ²

Stimulation Specifications:

Waveforms:	Interferential, Premodulated, Medium Frequency (Russian)
Output:	Interferential, 0-65 mA RMS, 1 kOhm load Premodulated, 0-50 mA RMS, 1 kOhm load Medium Frequency, 0-50 mA RMS, 1 kOhm load
Frequency:	Interferential, 4000 Hz, 0-250 Hz Premodulated, 4000 Hz, 1-250 Hz Medium Frequency, 2500 Hz, 50 Hz
Output Options:	Continuous, Vector Rotation, Surge, Reciprocation



Mettler Electronics Corp., 1333 S. Claudina St., Anaheim CA 92805 U.S.A.
(800) 854-9305 Department J Fax (714) 635-7539
Tel. (714) 533-2221 <http://www.mettlerelectronics.com>

