



Wave Goodbye To Dirt.

Now there's a way to safely wash away dirt and debris with ease. Let Cavitator® Ultrasonic Cleaners do the work for you.

Cavitator Cleaners generate ultrasonic waves that create millions of microscopic cleaning bubbles through the process of cavitation. The bubbles gently clean and polish any surface they come in contact with, penetrating every nook and cranny to dislodge trapped dirt and debris. Medical instruments, glassware, even jewelry are quickly restored to their



original lustrous finish, without expending any elbow grease.

The gentle ultrasonic cleaning action is ideal for delicate items, too.

Available in a variety of tank sizes, Cavitator Ultrasonic Cleaners are perfect for almost any size cleaning application. And they're simple to operate. Just fill with water, add a little Cavi-Clean detergent, and set the timer.

If you'd like to wave goodbye to dirt — and hard work — call 1-800-854-9305 for more information about the complete line of Cavitator Ultrasonic Cleaners.



METTLER **E**LECTRONICS®

Cavitator® Ultrasonic Cleaners

Cavitator Ultrasonic Cleaners are constructed from the highest quality stainless steel for extra-long life. Each model features an automatic timer for easy operation, plus solid state circuitry for exceptional reliability.

Model ME 18

Designed for high volume usage, like OR's, Central Supply or Labs, the Model ME 18 features a full 18 gallon capacity and is ideal for cleaning instruments or glassware prior to sterilization. Lockable, rolling rubber casters make it easy to move and store.




Model ME 5.5S

Our most popular size, the 5.5 gallon Model ME 5.5S is the perfect choice for medium duty applications. It can accommodate parts up to 21" in length, and is well-suited for almost any hospital or nursing home department, medical office, or outpatient surgery center.



The Smaller Cavitators


For a variety of lower volume applications, Cavitators are available in 2.1, 4.6 and 11 quart capacities. They fit neatly on any counter top and are perfect for countless medical, industrial and commercial applications.

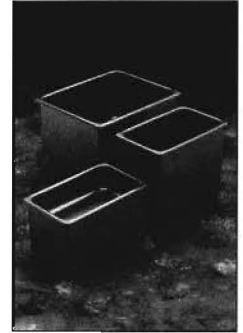


Specifications:

Model	ME 2.1	ME 4.6	ME 11	ME 5.5S	ME 18
Capacity	2.1 quarts (2 liters)	4.6 quarts (4.3 liters)	11 quarts (10.4 liters)	5.5 gallons (20.8 liters)	18 gallons (68 liters)
Outside Dimensions					
L x W x H					
inches	10 3/8 x 6 3/8 x 5 11/16	10 3/8 x 6 3/8 x 9 1/8	12 3/4 x 10 1/2 x 9	21 1/2 x 13 1/2 x 11 1/4	23 3/4 x 17 3/4 x 31
cm	26.4 x 16 x 14.4	26.4 x 16.2 x 23.2	32.4 x 26.7 x 22.9	54.6 x 34.3 x 28.6	60.3 x 45.1 x 78.7
Tank Size					
L x W x H					
inches	9 3/8 x 5 3/8 x 2 5/16	9 3/8 x 5 3/8 x 6	11 3/4 x 9 1/2 x 6	19 1/2 x 11 1/2 x 6	21 3/4 x 15 3/4 x 12
cm	23.8 x 13.7 x 5.9	23.8 x 13.7 x 15.2	29.8 x 24.1 x 15.2	49.5 x 29.2 x 15.2	55.2 x 40 x 30.5
Power Requirements	105/120 VAC 50/60 Hz	105/120 VAC 50/60 Hz	105/120 VAC 50/60 Hz	105/120 or 220 VAC 50/60 Hz	105/120 or 220 VAC 50/60 Hz
Power	85 watts	85 watts	200 watts	200 watts	800 watts
Nominal Main Frequency	67 KHz	67 KHz	67 KHz	67 KHz	24 KHz
Standard Accessories			Drain hose	Drain hose, Cleaning basket	Drain hose, Cleaning basket Fill hose, #1159

Accessories:

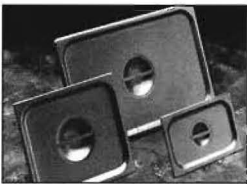

A variety of covers are available to help contain biologically-active contaminants.

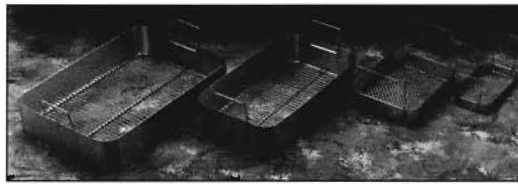

Baskets are available for each model to facilitate retrieval of cleaned items.



For optimum results, the Cavitator spray-rinse attachment is recommended to loosen dirt and debris prior to cleaning, and for a final rinse after cleaning.

Cavi-Clean is an all-purpose, biodegradable detergent specially formulated for use with Cavitator Ultrasonic Cleaners. It dissolves blood and inhibits rust formation, plus its neutral pH is gentle to the hands. Concentrated for economical use, dilute 1 part Cavi-Clean with 100 parts water for best results. For tougher jobs, Cavi-Clean Additive will boost the cleaning power of Cavi-Clean detergent.



METTLER ELECTRONICS®
corp.

1333 S. Claudina St., Anaheim, CA. 92805 (800) 854-9305 Department J
(714) 533-2221 FAX: (714) 635-7539